

PNGase F (Glycerol-free)



1-800-632-7799
info@neb.com
www.neb.com



P0705S 042151217121

P0705S




15,000 units Lot: **0421512** Exp: **12/17**
500,000 U/ml Store at **4°C** Do not freeze

Description: Peptide: N-Glycosidase F, also known as PNGase F, is an amidase which supplied glycerol free for optimal performance in HPLC intensive methods. PNGase F cleaves between the innermost GlcNAc and asparagine residues of high mannose, hybrid, and complex oligosaccharides from N-linked glycoproteins (1).

Source: PNGase F is purified from *Flavobacterium meningosepticum* (2).

Specificity:

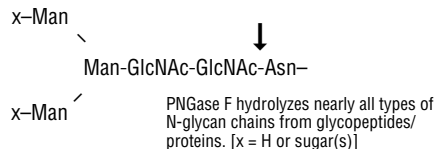

Applications:

- Removal of N-linked glycans from glycoproteins
- Preferred formulation for HPLC intensive methods

Supplied in: 50 mM NaCl, 20 mM Tris-HCl (pH 7.5 @ 25°C) and 5 mM Na₂EDTA.

Reagents Supplied with Enzyme:

10X Glycoprotein Denaturing Buffer:
(5% SDS, **0.4 M DTT**)

10X GlycoBuffer 2:
[0.5 M Sodium Phosphate (pH 7.5 @ 25°C)]

10% NP-40

Optimal incubation times and enzyme concentrations must be determined empirically for a particular substrate.

Reaction Conditions:

Typical reaction conditions are as follows:

1. Combine 1–20 µg of glycoprotein, 1 µl of 10X Glycoprotein Denaturing Buffer and H₂O (if necessary) to make a 10 µl total reaction volume.
2. Denature glycoprotein by heating reaction at 100°C for 10 minutes.
3. Make a total reaction volume of 20 µl by adding 2 µl 10X GlycoBuffer 2, 2 µl 10% NP-40, H₂O and 1–5 µl PNGaseF.
4. Incubate reaction at 37°C for 1 hour.

Note: Reactions may be scaled-up linearly to accommodate larger reaction volumes.

Unit Definition: One unit is defined as the amount of enzyme required to remove > 95% of the carbohydrate from 10 µg of denatured RNase B in 1 hour at 37°C in a total reaction volume of 10 µl (65 NEB units = 1 IUB milliunit).

Unit Definition Assay: 10 µg of RNase B are denatured with 1X Glycoprotein Denaturing Buffer at 100°C for 10 minutes. After the addition of NP-40 and GlycoBuffer 2, two-fold dilutions of PNGase F are added and the reaction mix is incubated for 1 hour at 37°C. Separation of reaction products are visualized by SDS-PAGE.

Quality Assurance: No contaminating exoglycosidase or Endoglycosidase F₁, F₂ or F₃ activity could be detected. No contaminating proteolytic activity could be detected.

Molecular Weight: 36,000 daltons.

Quality Controls

Glycosidase Assays: 5,000 units of PNGase F were incubated with 0.1 mM of fluorescently-labeled oligosaccharides and glycopeptides, in a 10 µl reaction for 20 hours at 37°C. The reaction products were analyzed by TLC for digestion of substrate.

(See other side)

CERTIFICATE OF ANALYSIS

PNGase F (Glycerol-free)



1-800-632-7799
info@neb.com
www.neb.com



P0705S 042151217121

P0705S




15,000 units Lot: **0421512** Exp: **12/17**
500,000 U/ml Store at **4°C** Do not freeze

Description: Peptide: N-Glycosidase F, also known as PNGase F, is an amidase which supplied glycerol free for optimal performance in HPLC intensive methods. PNGase F cleaves between the innermost GlcNAc and asparagine residues of high mannose, hybrid, and complex oligosaccharides from N-linked glycoproteins (1).

Source: PNGase F is purified from *Flavobacterium meningosepticum* (2).

Specificity:

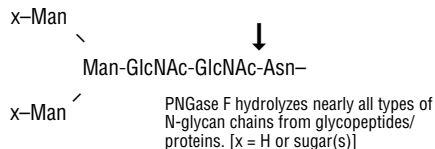

Applications:

- Removal of N-linked glycans from glycoproteins
- Preferred formulation for HPLC intensive methods

Supplied in: 50 mM NaCl, 20 mM Tris-HCl (pH 7.5 @ 25°C) and 5 mM Na₂EDTA.

Reagents Supplied with Enzyme:

10X Glycoprotein Denaturing Buffer:
(5% SDS, **0.4 M DTT**)

10X GlycoBuffer 2:
[0.5 M Sodium Phosphate (pH 7.5 @ 25°C)]

10% NP-40

Optimal incubation times and enzyme concentrations must be determined empirically for a particular substrate.

Reaction Conditions:

Typical reaction conditions are as follows:

1. Combine 1–20 µg of glycoprotein, 1 µl of 10X Glycoprotein Denaturing Buffer and H₂O (if necessary) to make a 10 µl total reaction volume.
2. Denature glycoprotein by heating reaction at 100°C for 10 minutes.
3. Make a total reaction volume of 20 µl by adding 2 µl 10X GlycoBuffer 2, 2 µl 10% NP-40, H₂O and 1–5 µl PNGaseF.
4. Incubate reaction at 37°C for 1 hour.

Note: Reactions may be scaled-up linearly to accommodate larger reaction volumes.

Unit Definition: One unit is defined as the amount of enzyme required to remove > 95% of the carbohydrate from 10 µg of denatured RNase B in 1 hour at 37°C in a total reaction volume of 10 µl (65 NEB units = 1 IUB milliunit).

Unit Definition Assay: 10 µg of RNase B are denatured with 1X Glycoprotein Denaturing Buffer at 100°C for 10 minutes. After the addition of NP-40 and GlycoBuffer 2, two-fold dilutions of PNGase F are added and the reaction mix is incubated for 1 hour at 37°C. Separation of reaction products are visualized by SDS-PAGE.

Quality Assurance: No contaminating exoglycosidase or Endoglycosidase F₁, F₂ or F₃ activity could be detected. No contaminating proteolytic activity could be detected.

Molecular Weight: 36,000 daltons.

Quality Controls

Glycosidase Assays: 5,000 units of PNGase F were incubated with 0.1 mM of fluorescently-labeled oligosaccharides and glycopeptides, in a 10 µl reaction for 20 hours at 37°C. The reaction products were analyzed by TLC for digestion of substrate.

(See other side)

CERTIFICATE OF ANALYSIS

No other glycosidase activities were detected (ND) with the following substrates:

β-N-Acetyl-glucosaminidase:

GlcNAcβ1-4GlcNAcβ1-4GlcNAc-AMC ND

α-Fucosidase:

Fucα1-2Galβ1-4Glc-AMC Galβ1-4
(Fucα1-3)GlcNAcβ1-3Galβ1-4Glc-AMC ND

β-Galactosidase:

Galβ1-3GlcNAcβ1-4Galβ1-4Glc-AMC ND

α-Galactosidase:

Galα1-3Galβ1-4Galα1-3Gal-AMC ND

α-Neuraminidase:

Neu5Acα2-3Galβ1-3GlcNAcβ1-3
Galβ1-4Glc-AMC ND

α-Mannosidase:

Manα1-3Manβ1-4GlcNAc-AMC
Manα1-6Manα1-6(Manα1-3)Man-AMC ND

Page 2 (P0705)

β-Glucosidase:

Glcβ1-4Glcβ1-4Glc-AMC ND

β-Xylosidase:

Xylβ1-4Xylβ1-4Xylβ1-4Xyl-AMC ND

β-Mannosidase:

Manβ1-4Manβ1-4Man-AMC ND

Endo F₁, F₂, H:

Dansylated invertase high mannose. ND

Endo F₂, F₃:

Dansylated fibrinogen biantennary. ND

Endoglycosidase F1 Assay: After incubation of 5,000 units of PNGase F with 20 pmol of 2-AA Man-5 fluorescent standard, for 20 hours at 37°C, no endoglycosidase F1 activity could be detected by LC/MS analysis with fluorescence detection.

Protease Assay: After incubation of 10,000 units of PNGase F with 0.2 nmol of a standardized mixture of proteins, for 20 hours at 37°C, no proteolytic activity could be detected by SDS-PAGE.

Physical Purity: Purified to > 95% homogeneity as determined by SDS-PAGE analysis using Coomassie Blue detection.

Heat Inactivation: 500 units of enzyme were inactivated by incubation at 75°C for 10 minutes.

Notes: Since PNGase F activity is inhibited by SDS, it is essential to have NP-40 present in the reaction mixture. Why this non-ionic detergent counteracts the SDS inhibition is unknown at present.

To deglycosylate a native glycoprotein, longer incubation time as well as more enzyme may be required.

PNGase F will not cleave N-linked glycans containing core α1-3 Fucose.

Previously supplied as a recombinant.

Repeated freeze thaw cycles degrade enzyme activity over time.

References:

1. Maley, F. et al. (1989) *Anal. Biochem.* 180, 195–204.
2. Plummer, T.H., Jr. and Tarentino, A.L. (1991) *Glycobiology* 1, 257–263.

Companion Product:

RNase B (NEB #P7817S)



NEW ENGLAND BIOLABS® is a registered trademark of New England Biolabs, Inc.

This product is intended for research purposes only. This product is not intended to be used for therapeutic or diagnostic purposes in humans or animals.

No other glycosidase activities were detected (ND) with the following substrates:

β-N-Acetyl-glucosaminidase:

GlcNAcβ1-4GlcNAcβ1-4GlcNAc-AMC ND

α-Fucosidase:

Fucα1-2Galβ1-4Glc-AMC Galβ1-4
(Fucα1-3)GlcNAcβ1-3Galβ1-4Glc-AMC ND

β-Galactosidase:

Galβ1-3GlcNAcβ1-4Galβ1-4Glc-AMC ND

α-Galactosidase:

Galα1-3Galβ1-4Galα1-3Gal-AMC ND

α-Neuraminidase:

Neu5Acα2-3Galβ1-3GlcNAcβ1-3
Galβ1-4Glc-AMC ND

α-Mannosidase:

Manα1-3Manβ1-4GlcNAc-AMC
Manα1-6Manα1-6(Manα1-3)Man-AMC ND

Page 2 (P0705)

β-Glucosidase:

Glcβ1-4Glcβ1-4Glc-AMC ND

β-Xylosidase:

Xylβ1-4Xylβ1-4Xylβ1-4Xyl-AMC ND

β-Mannosidase:

Manβ1-4Manβ1-4Man-AMC ND

Endo F₁, F₂, H:

Dansylated invertase high mannose. ND

Endo F₂, F₃:

Dansylated fibrinogen biantennary. ND

Endoglycosidase F1 Assay: After incubation of 5,000 units of PNGase F with 20 pmol of 2-AA Man-5 fluorescent standard, for 20 hours at 37°C, no endoglycosidase F1 activity could be detected by LC/MS analysis with fluorescence detection.

Protease Assay: After incubation of 10,000 units of PNGase F with 0.2 nmol of a standardized mixture of proteins, for 20 hours at 37°C, no proteolytic activity could be detected by SDS-PAGE.

Physical Purity: Purified to > 95% homogeneity as determined by SDS-PAGE analysis using Coomassie Blue detection.

Heat Inactivation: 500 units of enzyme were inactivated by incubation at 75°C for 10 minutes.

Notes: Since PNGase F activity is inhibited by SDS, it is essential to have NP-40 present in the reaction mixture. Why this non-ionic detergent counteracts the SDS inhibition is unknown at present.

To deglycosylate a native glycoprotein, longer incubation time as well as more enzyme may be required.

PNGase F will not cleave N-linked glycans containing core α1-3 Fucose.

Previously supplied as a recombinant.

Repeated freeze thaw cycles degrade enzyme activity over time.

References:

1. Maley, F. et al. (1989) *Anal. Biochem.* 180, 195–204.
2. Plummer, T.H., Jr. and Tarentino, A.L. (1991) *Glycobiology* 1, 257–263.

Companion Product:

RNase B (NEB #P7817S)



NEW ENGLAND BIOLABS® is a registered trademark of New England Biolabs, Inc.

This product is intended for research purposes only. This product is not intended to be used for therapeutic or diagnostic purposes in humans or animals.